www.artemismask.com

LED MASK User' s Manual

Product Name : miin-iMask Artemis, Artemis II

產品名 : miin-iMask Artemis, Artemis II

Model No : Artemis-F

型號名 :Artemis-F

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Please read the user manual carefully before using the device

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Safety Information

1. Caution for the device

- Do not use the device if you find smoke or strange smell from it and power it off immediately.
- Use the power adapter and USB cable provided.
- Do not disassemble, repair, modify.
- Do not hold the power cord with wet hands.
- Remove the power plug during thunderstorm.
- Remove the power plug when carrying the device.
- Use the device in the place without electromagnetic interference.
- Do not drop or shock the device
- Do not use the device if it is dropped in the water and contact the customer center.
- Do not prick the unit gap with sharp metal objects such as keys, hairpins, necklaces, or nails.
- Replace the damaged device or accessories such as power adapter with genuine products.
- The warranty service is not available if you disassemble the device.
- Do not use the power cord with damaged insulator.
- Do not scratch or bend the cord severely.
- Do not copy the design, software without permission.

2. Caution for injury

- Do not use the device for purposes other than those set out in this guide.
- The device is recommended to be used once a day.
- Do not make direct eye contact with LED light while device is in operation. Excessive exposure to LED light can cause eye damage.
- Use the device with eye closed or eye patch if your eyes are sensitive to the light.
- Do not use the device if you have skin disease or skin injury .
- Do not use this device if you are currently undergoing or underwent in the last 2 months, skin related procedures, cosmetic surgery or procedures for hair loss. Please only use it after consultation with a specialist.
- Use the device after consultation with a dermatologist if you have allergic or sensitive skin.
- The device may cause skin irritation such as redness. In this case, use the device after consultation with a dermatologist.
- use the device after consultation with a specialist if any of the following applies:
 - Cardiovascular disease, Endocrine abnormality, Sensory nervous system abnormality, Skin disease
 - Ophthalmology and surgery (LASIK, LASEK, etc.), Infectious disease, A person with epileptic symptoms
 - A person taking hormone imbalance medication
- Do not leave in the car under the hot sun or near heat sources. There is a risk of burns.
- Do not store or use in hot and humid place.
- Do not allow children to play with the product.
- Pregnant, menstrual, lactating women are not recommended for use. Unstable hormonal balance can lead to severe skin problems.
- Do not use if you have photosensitive diseases such as sunlight allergy.
- Do not use while driving or walking.
- Do not use if you have implanted an electronic device such as a pacemaker or have a medical electronic device such as a medical monitoring device.



3. Caution for battery

- The damaged battery has the risk of explosion. Do not shock or disassemble the battery
- Remove power cord from the device if you do not use it for a long time.
- The battery is not charged when you purchase the device. Charge it before using it.


4. Product Storage and maintenance

- Do not store the device in the high place. You may be injured off.
- Keep the device away from heated place. It may cause fire.
- Do not drop the device or the charging adapter into a drain or sink or place it in a place where it may fall into water or other liquids.
- Keep the cable away from heated surface.
- Do not place volatile liquids, such as pesticides, near them
- Wiping the unit surface hard can damage the unit.
- Do not place plastic products or rubber on the unit for a long time.
- Do not put the unit in dishwasher, oven, washing machine.
- To clean the unit, use a soft, dry cloth. If the surface is too dirty, use a soft cloth with water or a wet tissue. Do not use strong chemicals such as alcohol, benzene or thinner as they may deteriorate the surface of the unit.
- Do not leave the device in the car under the sun or near the heat place.
- Keep the device away from the hot and humid place.

Parts / Unit Description

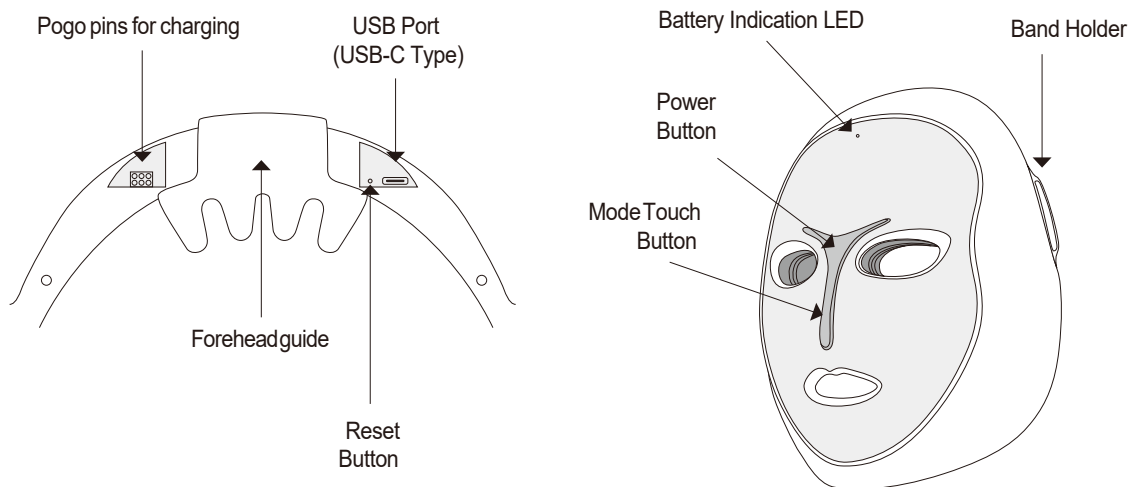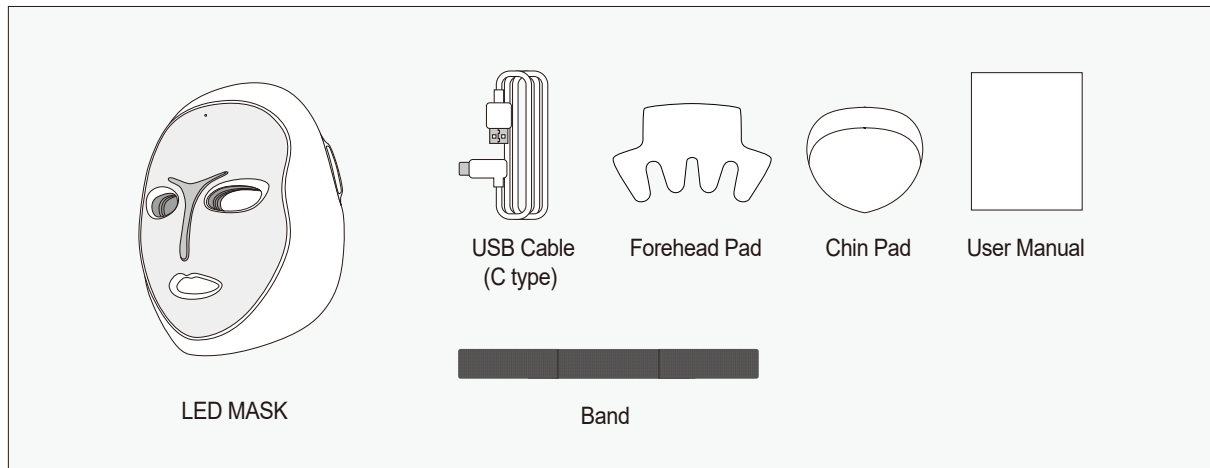Charging battery

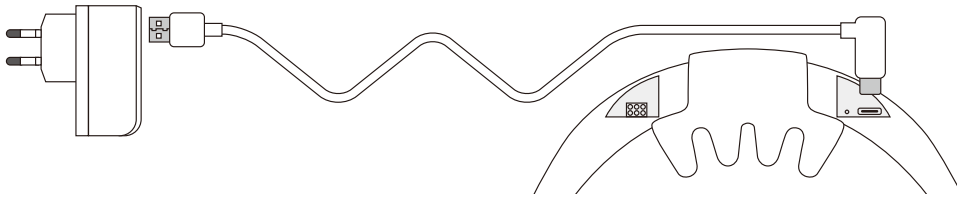
1. Connect USB cable to battery charging adapter.
2. Plug the adapter into the electric outlet.
3. Connect USB cable to the mask.
4. LED color changes, depending on the remaining battery.
 - The battery charging time is subject to the temperature, remaining battery or charging condition.
 - The battery charging time is around 3.5 hours.
 - The device can be used while charging the battery.
 - The battery should be charged if LED color is red.

Battery Status	LED Color
Fully charged	Blue
High	Blue Flickering
Middle	Green Flickering
Low	Red Flickering
Need to charge	Red Flickering


Before using the device

Use the foundation cosmetics after washing the face

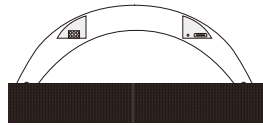Adjust the band to fit your face to the mask

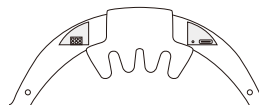
Use the rubber pad If your face is small

Wear the mask from top to bottom


Using the device

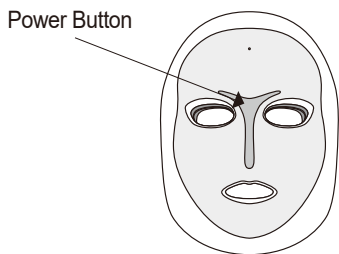Power Button

Power On

Long press the power button
Press for more than 1 second to power on

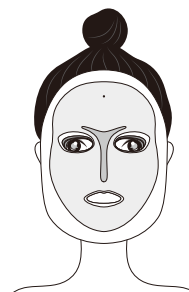
Skincare starts

3 seconds after power on
Start management automatically

1. Power on by pressing the power button for 1 second. (RGB LED turns on in order.)
2. Skincare mode automatically starts after 3 seconds. (Default setting : Normal mode (red color))
3. Skincare mode stops after 10 minutes.(Default settings: 10 minutes)
4. The power is off 3 seconds after the skincare mode stops.
5. The power is off by pressing the power button for 1 second in PAUSE mode or SKINCARE mode.

* The time and intensity is adjustable in the app.


Mode Touch button
(Touch briefly)

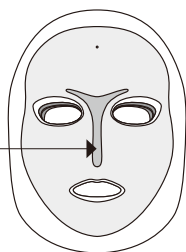Changing Normal Mode

If touch the touch button briefly during operation:
change to normal mode
(Red → Green → Blue → Yellow → Purple → Sky Blue)

Mode Touch button
(Touch long)

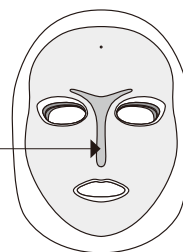Changing Auto Mode


When touching touch button for longer than 1 second
during operation: Auto mode change
(Auto 1 → Auto 2 → Auto 3 → Auto 4 → Auto 5 → Auto 6 → Auto 7)

1. Changing Normal Mode
 - Change the normal mode by pressing the touch button briefly
 - Red → Green → Blue → Yellow → Purple → Sky Blue
2. Changing Auto Mode
 - Change the normal mode by pressing the touch button over 1 second.
 - Whenever touch the button over 1 second.
Auto Mode 1 → Auto Mode 2 → Auto Mode 3 → Auto Mode 4 → Auto Mode 5 → Auto Mode 6 → Auto Mode 7 → Auto Mode 1

* Other care modes are convenient to select using the app.

Choose a mode

Mode Button
(Touch briefly)

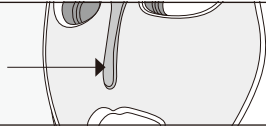The mode changes whenever touching the button.
(Red → Green → Blue → Yellow → Purple → Sky Blue)

Normal Mode 1
(RED)

Elasticity,
Antiaging


Normal Mode 2
(GREEN)

Skin Tone-up

Normal Mode 3
(BLUE)

Oily Skin

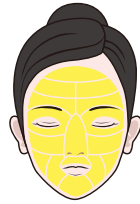Normal Mode 4
(YELLOW)

Reducing
Redness,
Whitening

Normal Mode 5
(VIOLET)


Deep cleansing,
Healing after the
treatment

Normal Mode 6
(SKY BLUE)

Moisturizing

Mode Button
(Touch long)

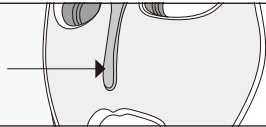The mode changes whenever touching the button for 1 second.
Auto mode 1 → Auto mode 2 ••• Auto mode 7 → Auto mode 1

Auto Mode 1
(Elasticity Care)

Elasticity,
regeneration care

Auto Mode 2
(Trouble Care)

Acne immediately
after procedure

Auto Mode 3
(Whitening Care)

Whitening, flushing care


Auto Mode 4
(Skin Soothing, Moisturizing)

Sensitive skin,
water balance


Auto Mode 5
(Homme Care)

Aftershave,
elasticity care

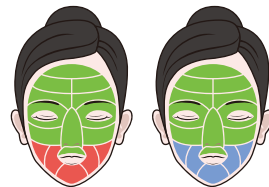Auto Mode 6
(V-line Care)

Chin line intensive care

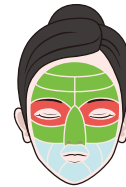Auto Mode 7
(Fine dust care)

Fine dust, O2 care

Installing app

Android	<ul style="list-style-type: none"> •Play Store → Search “iMask” → Install •Register the device in the app. 	 
iOS	<ul style="list-style-type: none"> •App Store → Search “iMask” → Install •Register the device in the app . 	 
Bluetooth Connection	<ul style="list-style-type: none"> •Bluetooth is activated when the power is on . •Select “miin-iMask-xxxx” in the Bluetooth list of your smartph one to connect 	
Smartphone Minimum Specifications for App.	<ul style="list-style-type: none"> •Android Phone : OS – Android 6.0 or higher / Display – 1080 x 1920 or higher •iPhone : OS – iOS 11.0 or higher / Display – 750 x 1334 or higher 	

Troubleshooting

Trouble	Solution
No power	<ul style="list-style-type: none"> • The battery may be empty. Charge the battery. • The device can be used while charging the battery. • Check the right connection between the product and the electric outlet • Check if the charger is working. • Check if the voltage is correct. • Contact the customer center if the battery runs out quickly.
The device is not working correctly.	<ul style="list-style-type: none"> • Press the power button over 1 second. ** The power button is not the touch one. • Do not hold the touch button but touch it. (Short touch : Changing the normal mode, Long touch : Changing the auto mode) • Use the device after charging the battery. • Check if the battery indication LED is on
Battery indication LED is not on	<ul style="list-style-type: none"> • Check if USB cable is connected to the charger correctly. • Check if the charger is connected to the electric outlet. • Use other USB cable • Contact the customer center if the problem is not solved.
The device is heating	<ul style="list-style-type: none"> • The device can be heating while it is charged or used. It is normal and do not affect the device life.

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specification

Product Name	miin-iMask	LED	RGB LED, IR LED
Model No.	Artemis-F	Operation Time	Default Settings : 10 minutes (Adjust 1~30 minutes)
Weight	429 g	Continuous Operation Time	2 hours
size	205 (W) X 260 (D) X 135 (H) mm	Battery full charging time	3.5 hours
Power Adapter	Input : 100~240V, 50/60Hz Output : 5V, 2A	Operation temperature / Humidity	Temperature : 0°C~35°C, Humidity : 0%~90%
Battery	3.7V, Li-ion Rechargeable Battery, 3000mAh / Artemis II 3700 mAh	Recommended Care Time	10 minutes / Daily

■ Responsibility for paid service (customer burden)

When applying for service, the following cases will be charged for free even during the free period.

- If a malfunction occurs due to foreign substances in the product (water, beverages, coffee, toys, etc.)
- In case of product damage or malfunction due to dropping or impact
- In case of appearance damage or deformation caused by thinner, benzene or other organic solvent
- When not using the original product (adapter, USB cable, etc.)
- In case of product failure due to incorrectly applied voltage
- If the customer disassembles the product arbitrarily and the accessory is lost or damaged.

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