# iNeck - User's Manual (CHN)

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

# iNeck - User's Manual (ENG)

## ■ Warnings

⚠ If you do not follow the instructions, it may cause equipment damage and personal injury.

- Do NOT put the device in the heater or microwave range. It may cause device damage or explosion.
- Do NOT use or charge when it gets wet. It may become the cause of a fire, leaks, excessive heat or rupture.
- NEVER turn it ON when it gets wet. It may cause malfunction.
- Do NOT disassemble or give the device a shock. It may be damaged if you drop it on the hard floor.
- Keep your face or hands from the device when it is broken. It may cause your skin injury.
- Please put in the power plug tightly. Do NOT bend, twist, heat or cut the power cord. It can be broken and cause a fire or electronic shock.
- Do not look into the LED while the device is in operation. Your eyesight may be impaired.
- Do not remove the product while the device is in operation.
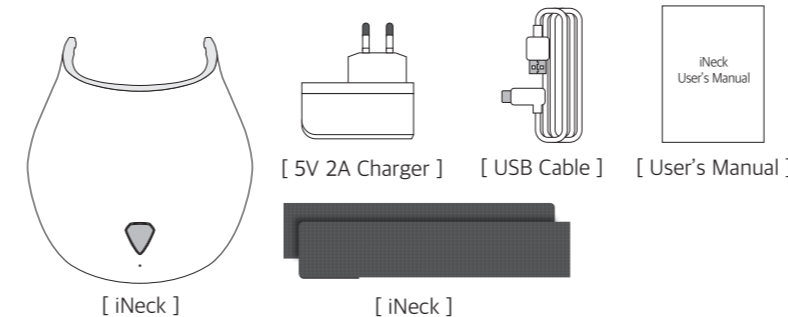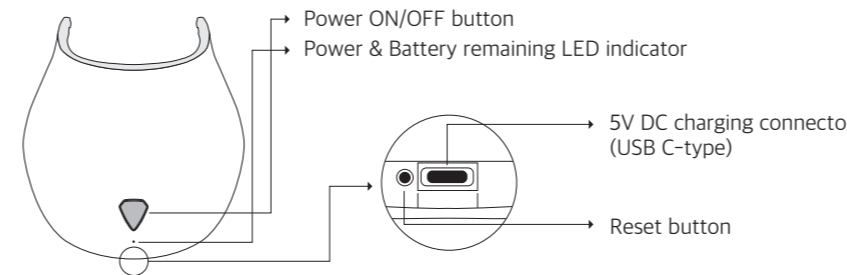
## ■ Power ON

- Power on with the beep sound to press the power button for 1 second.
- LED indicator is on.
- The skincare management mode automatically starts in 3 seconds after the power is on.

## ■ Automatic Power OFF

- The power automatically turns off in 3 seconds after the skincare is finished

## ■ Power OFF

- Power off with the beep sound to press the power button for 1.5 seconds.

## ■ Charging the battery

- Connect the device to 5V DC adapter with USB cable.
- The full charging time is 3 hours

## ■ Mode Changing

- The skin management mode automatically starts in 3 seconds when the power is on.
- Change the mode by pressing the power button shortly during the skincare mode.
  - : Mode 1 → Mode 2 → Mode 3 → Mode 4 → Mode 5 → Mode 1 repeat
  - Mode 1 : Red
  - Mode 2 : Green
  - Mode 3 : Blue
  - Mode 4 : Red + Blue
  - Mode 5 : Blue + Red

# iNeck - User's Manual (ENG)

## ■ Before using the device

\*Use the foundation cosmetics after washing the face

Please wash your neck. Wear the iNeck. Adjust the band to fit your neck to the iNeck

## ■ Install App. & Authentication Key Registration

- Download the app from Apps store of Android & iOS smart phone.
- Play Store / App Store → Search "iMask" → Install

- Use the app after the product registration
- After running the app, you will be prompted to enter your product activation key number
- Enter "1 2 3 4 5 6 7 8 9 0" (10 digits) in order.
- (Caution) Smartphone Minimum Specifications for App.
 


|                                 |                                |
|---------------------------------|--------------------------------|
| [Android Phone]                 | [iPhone]                       |
| OS - Android 6.0 or higher      | OS - iOS 11.0 or higher        |
| Display - 1080 x 1920 or higher | Display - 750 x 1334 or higher |

## ■ Bluetooth connection

- Bluetooth is activated when the power turns on.
- On your smartphone, "miin-iMask Neck" appears in the bluetooth list as follows and it will be automatically connected if you select it

## ■ Battery indicator

| Step | State    | Working      | Charging       |
|------|----------|--------------|----------------|
| 1    | Very low | Blinking Red | Blinking Red   |
| 2    | Low      | Red          | Blinking Red   |
| 3    | Niddle   | Green        | Blinking Green |
| 4    | High     | Blue         | Blinking Blue  |
| 5    | Full     | Blue         | Blue           |
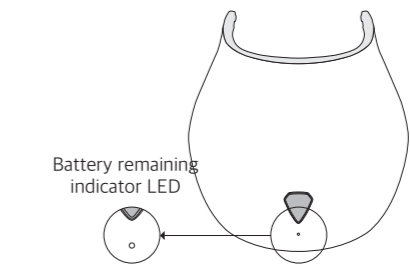
## ■ Specifications

|              |  |
|--------------|--|
| Product Name | miin-iMask   |
| Model Name   | Artemis-N  |
| Weight       | 287g   |
| Size         | 205(W) X 245(D) X 105(H)mm                                     |
| Components   | iMask-Neck, Neck Band, User's Manual, DC 5V charger, USB Cable |
| Power supply | 3.7V, Li-ion Rechargeable Battery, 3000mAh                     |

## ■ Warranty

1. This product is manufactured under strict product quality compliance and examination process of Miinsystems co., ltd.
2. Repair is free of charge only when malfunction happens under users' ordinary use(Within Warranty Period)
3. Repair needs to be paid, if the type of breakdown is out of Free A/S service.

|                    |            |                 |                          |
|--------------------|------------|-----------------|--------------------------|
| Product Name       | miin-iMask | Warranty Period | from the purchase 1 year |
| Model Name         | Artemis-N  |                 |                          |
| Serial No.         |            |                 |                          |
| Date of Sales      |            |                 |                          |
| Customer's Address |            |                 |                          |
| Customer's Name    |            | tel             |                          |

- Guide for Free repair In case the product is broken within 12 months from the purchase, free repair would be provided when it has been used in ordinary use environment.
- Guide for Paid repair Repair needs to be paid in the following cases.
  - Malfunction caused by user's mistake and mishandling failure.
  - Out of warranty period.
  - Malfunction caused by customer's intention or mistake(carelessness).
  - Malfunction and damage caused by other service engineer's repair
- Product consultation and A / S guide : www.miinsystems.com / +82-31-385-9969